

BioSizer-100™ DNA Ladder Plus

Catalogue Number	Description	Qty
BG-M100	BioSizer-100 DNA Ladder Plus	50ug

Bio100 DNA Ladder contains 13 linear double strand DNA bands of 100, 200, 300, 400, 500, 600, 700, 800, 900, 1000, 1500, 2000 and 3000 bases pairs. The intensity of 500, 1500, 2000, 3000bp bands are stronger than others, for instant band identification.

This marker can also be used as reference for DNA quantitative analysis.

DNA quantity of each band per 5µl loading are as follows:

- 💧 50ng for 100, 200, 300, 400, 600, 700, 800 and 900bp bands
- 💧 130ng for 500bp band
- 💧 100ng for 1000bp band
- 💧 200ng for 1500, 2000 and 3000bp bands


ANALYSIS

The marker is evaluated for resolution, intensity, and background on 2% agarose Gel in 0.5X TAE buffer.


STORAGE CONDITION


-20° C, stable for one year


Bio/100 DNA Ladder

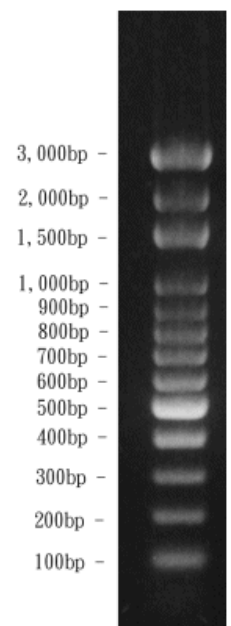


2% Agarose gel in 0.5X TAE Buffer

 BioGene Limited – 6 The Business Centre, Harvard Way, Kimbolton, Cambridgeshire, PE28 0NJ, United Kingdom.

 +44 (0)3333 222678

 biogene@biogene.com

 +44 (0)3333 222567

 www.biogene.com

Registered No. 2800327

VAT No. GB576854584



5282-QMS-001
ISO 9001



0026